

Import of Referral Lab Results to PowerChart

Requirements:

Locate accession in Cerner ORV


Locate results from reference lab web page or Quantum

Ability to apply FIN on Adobe PDF files


1. Locate results from referral lab web page or Quantum
 - a. Save results as PDF to secure folder (as designated by site management staff)
 - b. Apply "FIN#" (ex. FIN# 12345678) to uppermost right hand corner using adobe writer.
 - c. Save updated file to folder.



2. Open Cerner ORV and locate Accession to be result to complete.
3. Open Result Entry and apply "See Report" in the entry
 - a. Verify Result entry
 - b. Refresh Cerner to ensure accession is "Complete" status in ORV
4. Open PowerChart using the Patient FIN #
 - a. Open Clinical Notes folder

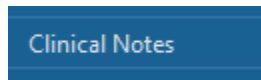

- i. Double-click on the scanner image to open the "add document" folder

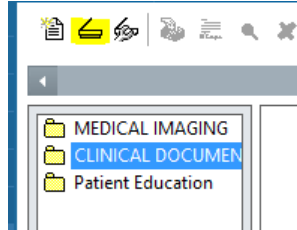


- c. Using the Dropdown – select the appropriate report "Type"



*Type: [Dropdown menu]

- e. For this example- we are using ".Lab Rprt"

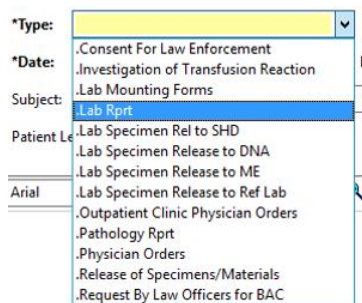

- g. Apply the correct date/time to match the accession information in ORV

Import of Referral Lab Results to PowerChart

*Date: 11/17/2021 1445 MST

- i.
- j. Apply the Primary name from Cerner as the "Subject"
 - i. Do not use the Cerner mnemonic as the subject name
 - ii. Apply the specimen type after the test name

Subject: Arbovirus Antibody Panel, IFA; CSF

- 1.
- k. Click on "Import" at the bottom of the screen

Paste File Sign Import Scan Save Save & Close Cancel

- i.
- ii. Select the report from saved folder (step 1)

File name: 43077734.pdf All Files (*.*)
Open Cancel

- 1.
- iii. The report will populate in the previously blank area of the screen

Types: Lab Rpt Patients: Spine Joints
Date: 10/30/2021 10:13 AM Status: In Progress
Subject: Arbovirus Antibody Panel, IFA, CSF
Patient Level Documents: ☐

SEND TO: BANNER-UNIV MED CNTR S CAMPUS
2800 E. AJO WAY
TUCSON, AZ 85713

Quest Diagnostics Nichols Institute™

FIN# 43077734
Quest Diagnostics, Incorporated
33608 Ortega Hwy., San Juan Capistrano, CA
CLIENT SERVICES - (800) 553-5445
Director: Nina Maramba, M.D., Ph.D.

1. ALIGN HERE

2. Verify accuracy of Patient to name on report

3. Verify all pages are appropriate to this order only

4. Verify all pages are for this patient.


5. Report with multiple pages with open to last page of the file imported

Page 1 of 2

- 6.
- iv. Once verification of correct patient and report; select Sign

Paste File Sign Import

- v.
- l. Results will load to "Clinical Documents" area; "Laboratory Reports"

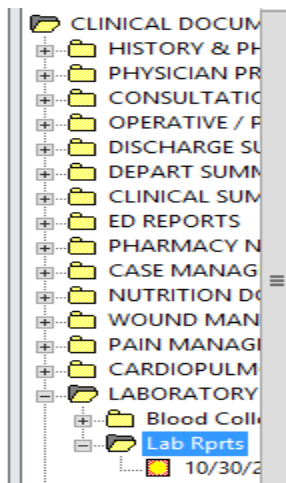

i.

Import of Referral Lab Results to PowerChart

- m. Verify results are posted to the correct patient.
 - i. Patient name in PowerChart vs report
 - ii. FIN on report
 - iii. All pages are legible and for the correct patient
 - iv. One Report per import